

Anti-Mouse CD279/PD1-159Tb

Catalog: 3159023B Clone: J43
Package Size: 100 tests Isotype: IgG


Storage: Store product at 4°C. Do not freeze. Formulation: Antibody stabilizer with 0.05% Sodium Azide

Reactivity: Mouse,

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 μ l for up to 3 X 10^6 live cells in 100 μ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

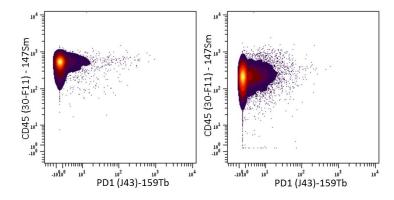

Mouse splenocytes were incubated for 3 days in media alone (left) or with Concanavalin A (right) and then stained with 147Sm-anti-CD45 (30-F11) and 159Tb-anti-PD1 (J43). T and B cells are displayed in the analysis.

Description

PD-1 (programmed cell death protein 1), also known as CD279, is a cell surface membrane protein member of the Ig superfamily. PD-1 contains the immunoreceptor tyrosine-based inhibitory motif (ITIM) and plays a key role in peripheral tolerance and autoimmune disease in mice. PD-1 is expressed mainly on activated T and B lymphocytes. The RMP1-30 antibody reacts with mouse PD-1.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly Multiparametric Analysis by Mass Cytometry. J Immunol Methods 361 (1-2):1-20, 2010.

$\label{prop:comport} \textbf{For technical support visit fluidigm.com/support}$

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies supplied by eBioscience.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2015 Fluidigm Corporation. All rights reserved. 07/2015