

# **Anti-Human Integrin-168Er**

Catalog: 3168013B Clone: Y9A2
Package Size: 100 tests Isotype: IgG1


Storage: Store product at 4°C. Do not freeze. Formulation: Antibody stabilizer with 0.05% Sodium Azide

Reactivity: Human, Guinea Pig,

## **Technical Information**

**Validation:** Each lot of conjugated antibody is quality control tested by  $\mathsf{CyTOF}^{\textcircled{R}}$  analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

**Recommended Usage:** The suggested use is 1  $\mu$ l for up to 3 X  $10^6$  live cells in 100  $\mu$ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

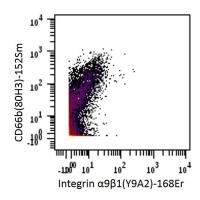

Human whole blood stained with 168Er-anti-Integrin  $\alpha 9\beta 1(Y9A2)$  and 152Sm-anti-CD66b(80H3). Total viable cells are displayed in the analysis.

## **Description**

The Y9A2 monoclonal antibody binds to human Integrin  $\alpha9\beta1$ , a transmembrane protein belonging to the integrin family. Multiple ligands bind to Integrin  $\alpha9\beta1$ , including the inducible endothelial counter-receptor, vascular cell adhesion molecule-1 (VCAM-1), the extracellular matrix proteins tenascin C and osteopontin, coagulation factor XIII, and von Willebrand factor. Integrin  $\alpha9\beta1$  is expressed on a variety of cell types including airway epithelium, different types of muscle cells, and hepatocytes. In blood, integrin  $\alpha9\beta1$  is highly expressed on human neutrophils and at low level on monocytes. Recent reports have found Integrin  $\alpha9\beta1$  expression on hematopoietic stem and progenitor cells, where it is thought to mediate stem cells adhesion, proliferation, and differentiation.

### References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly Multiparametric Analysis by Mass Cytometry. J Immunol Methods 361 (1-2):1-20, 2010.

### For technical support visit fluidigm.com/support

EMAIL

 $\label{lem:composition} United States support.northamerica@fluidigm.com \mid Europe support.europe@fluidigm.com \\ Asia techsupportasia@fluidigm.com \mid Latin America techsupporlatam@fluidigm.com \\ All other countries techsupport@fluidigm.com \\ \end{tabular}$ 

PHONE
United States (toll-free) +1 866 358 4354 | Europe +33 1 60 92 42 40
Japan +81 3 3662 2150 | China (excluding Hong Kong) +86 21 3255 8368
All other countries +1 650 266 6100

#### For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, Cell-ID, CyTOF, and Maxpar are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. All other trademarks are the sole property of their respective owners. © 2015 Fluidigm Corporation. All rights reserved. 06/2015