

Maxpar Antibody Conjugation Service for Mass Cytometry

Expand and accelerate your research with isotope-conjugated antibodies made by experts



Mass cytometry is the most comprehensive way to profile cell phenotype and function, enabling simultaneous detection of more than 40 cellular targets in every precious sample.


Whether you're building a 15- or 40-parameter mass cytometry panel, the extensive Fluidigm catalog of pre-conjugated antibodies and antibody panel kits will get you there.

We understand your panel design needs. If you have a proprietary antibody to add to your panel or require customized design, our Maxpar® Antibody Conjugation Service is committed to your success.

Select the desired isotope from our 37 high-purity metal tags, decide which confirmation method you prefer—biochemical analysis or more extensive biological functional testing—and place your order.

After our expert scientists perform the conjugation and analysis, we'll send your antibody conjugate along with a detailed analytic report.

Custom Conjugation Workflow

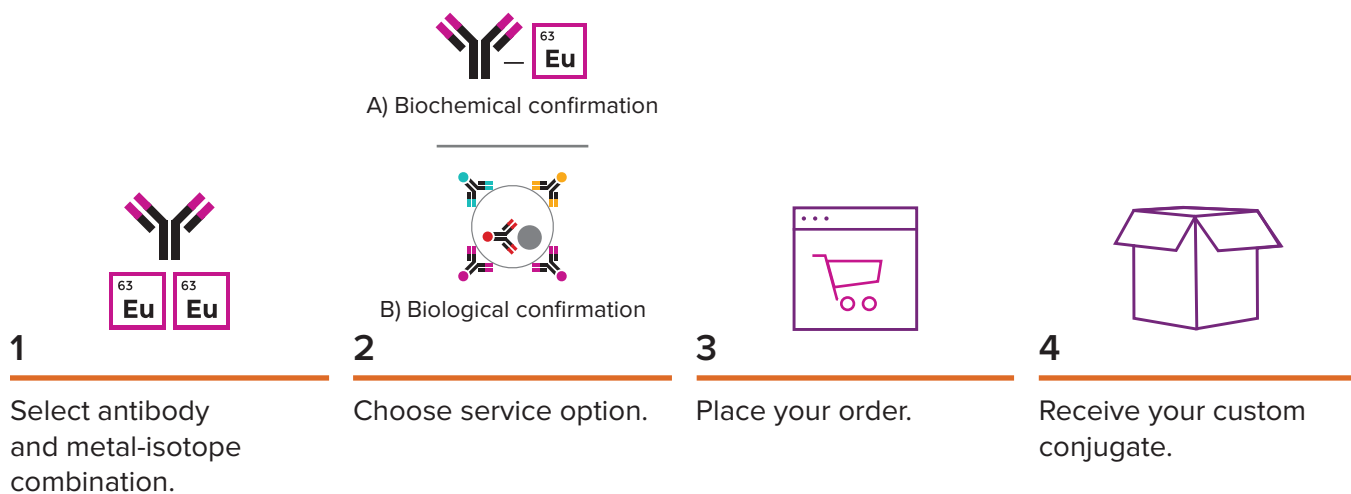

Highlights

Flexible—Extend your mass cytometry panels by adding metal-labeled targets not yet available in the Maxpar catalog.

Efficient—Save time and resources by letting the experts do the labeling and quality control for you.

Dependable—Conjugation service performed by our professional technicians ensures product with consistently high quality and reliable performance.

Custom Conjugation Service Specifications and Ordering Information

	Antibody Conjugation Service with Biological Confirmation	Antibody Conjugation Service with Biochemical Confirmation	
Biochemical analysis provided	Yes	Yes	
Biological confirmation performed	Yes (Available for Helios™/CyTOF only)*	No	
Staining profile and titer data	Yes	No	
Timeline	6 to 8 weeks from receipt of antibody	4 to 6 weeks from receipt of antibody	
Yield			
Standard scale Input 200 µg Purified, protein-free buffer	Cat. No. 3999999-1 Yield: 100 tests (1–3 x 10 ⁶ cells)	Cat. No. 3999999-2 Yield: 40 µg minimum	
Bulk scale Variable input mass of antibody Purified, protein-free buffer	Cat. No. 3999999-5 Extra 100-test vial Additional units of 100 tests each	Cat. No. 3999999-7 Yield: 200 µg minimum	Cat. No. 3999999-8 Yield: 600 µg minimum

*Not available for antibodies used on the Hyperion™ Imaging System

Ligands and buffer suitable for custom conjugation—IgG monoclonal and polyclonal antibodies are routinely conjugated. The feasibility of conjugating other antibody isotypes and proteins is considered on a case-by-case basis. To facilitate labeling, ligands must be provided purified, unconjugated and in a protein-free buffer at minimum concentration of 0.5 mg/mL.

Source of antibodies—Antibodies to be custom-conjugated may be either provided by the customer or purchased from a commercial source. Where Fluidigm purchases an antibody on behalf of a customer, the cost is quoted and billed as a separate line item, Cat. No. 3999999-3, Antibody Material Handling Charge.

Quality control—All custom conjugations are subject to rigorous biochemical confirmation, including inductively coupled plasma mass spectrometry, determination of concentration and yield and assessment of ratio of attachment of metal atoms to antibody. Optionally available biological confirmation service involves empirical testing by CyTOF® on a cell line, with optimal titer determined and the conjugate formulated for delivery in a 100-test vial.

To learn more about Maxpar Custom Conjugation Services contact your local Field Application Specialist or Sales Representative.

fluidigm.com/customconjugationservice

Corporate Headquarters

7000 Shoreline Court, Suite 100
South San Francisco, CA 94080 USA

Toll-free: 866 359 4354 in the US and Canada

Fax: 650 871 7152

For Research Use Only. Not for use in diagnostic procedures.

Information in this publication is subject to change without notice. **Patent and license information:** fluidigm.com/legalnotices. Fluidigm, the Fluidigm logo, CyTOF, Helios, Hyperion and Maxpar are trademarks and/or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2018 Fluidigm Corporation. All rights reserved. 09/2018

Sales

North America | +1 650 266 6170 | info-us@fluidigm.com

Europe/EMEA | +33 1 60 92 42 40 | info-europe@fluidigm.com

China (excluding Hong Kong) | +86 21 3255 8368
info-china@fluidigm.com

Japan | +81 3 3622 2150 | info-japan@fluidigm.com

All other Asian countries | +1 650 266 6170 | info-asia@fluidigm.com

Latin America | +1 650 266 6170 | info-latinamerica@fluidigm.com