

Anti-Human CD184/CXCR4-175Lu

Catalog: 3175001B

Package size: 100 tests

Storage: Store at 4 °C. Do not freeze.

Cross-reactivity: Cynomolgus Monkey, Rhesus, Chimpanzee, Olive Baboon, Sooty Mangabey, African Green

Clone: 12G5


Isotype: Mouse IgG2a

Formulation: Antibody stabilizer with 0.05% sodium azide

Technical Information

Validation: Each lot of conjugated antibody is quality control-tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended usage: The suggested use is 1 μ L for up to 3×10^6 live cells in 100 μ L. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

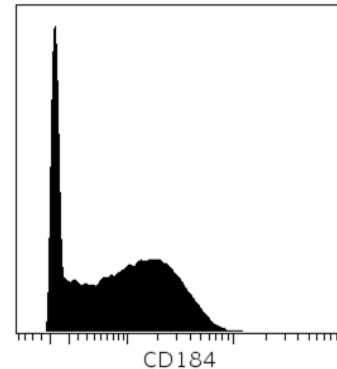

Human PBMC stained with 175Lu-anti-CD184/CXCR4 (12G5).

Description

CD184, also known as CXCR4, is an alpha-chemokine receptor specific for stromal-derived factor-1 (SDF-1). This surface protein has seven transmembrane regions. CXCR4 is one of the most important chemokine receptors that act as co-receptors for HIV entry into CD4+ T cells. It is expressed widely in blood and several other tissues.

References

Bandura, D.R. et al. "Mass cytometry: technique for real time single cell multitarget immunoassay based on inductively coupled plasma time-of-flight mass spectrometry." *Analytical Chemistry* 81 (2009): 6,813–22.

Ornatsky, O.I. et al. "Highly multiparametric analysis by mass cytometry." *Journal of Immunological Methods* 361 (2010): 1–20.

Bendall, S. C., et al. Single-Cell Mass Cytometry of Differential Immune and Drug Responses Across a Human Hematopoietic Continuum. *Science* (2011): 687–96.

For technical support visit <http://techsupport.fluidigm.com>. | For general support visit www.fluidigm.com/support.

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend[®] and licensees thereof. Information in this publication is subject to change without notice. **Safety data sheet information:** www.fluidigm.com/sds. **Patent and license information:** www.fluidigm.com/legalnotices. **EU's WEEE directive information:** www.fluidigm.com/compliance. Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. All other trademarks are the sole property of their respective owners. © 2017 Fluidigm Corporation. All rights reserved. 11/2017