

# Anti-Human CD223/LAG-3-153Eu

## Pathologist-Verified Clone for Imaging Mass Cytometry™

Catalog: 3153028D

Package size and concentration: 25 µg, 0.5 mg/mL

Storage: Store at 4 °C. Do not freeze.

Reactivity: Human

Clone: D2G40

Isotype: Rabbit IgG

Formulation: Antibody stabilizer with 0.05% sodium azide

Application: IMC-Paraffin

## Technical Information


**Application:** The metal-tagged antibody is designed and formulated for the application of Imaging Mass Cytometry (IMC™) using the Fluidigm Hyperion™ Imaging System on formalin-fixed, paraffin-embedded (FFPE) tissue sections.

**Quality control:** Each lot of conjugated antibody is quality control-tested by Imaging Mass Cytometry on tissue sections.

**Recommended concentration:** For optimal performance it is recommended that the antibody be titrated for the desired application. Suggested initial dilution range:  
IMC-Paraffin: 1:25 to 1:100

## Description

Lymphocyte activation gene 3 (LAG-3), also known as CD223, is a 70 kDa surface glycoprotein belonging to the Ig superfamily. LAG-3 shares homology with CD4, and LAG-3 binding to MHC class II is thought to contribute to the negative regulation of T cell activation. LAG-3 is expressed on T cells, including regulatory T cells and NK cells. Recently, it has been suggested that co-expression of CD49b and LAG-3 identifies human and mouse type 1 regulatory T (Tr1) cells. The D2G40™ monoclonal antibody recognizes endogenous levels of total LAG-3 protein.

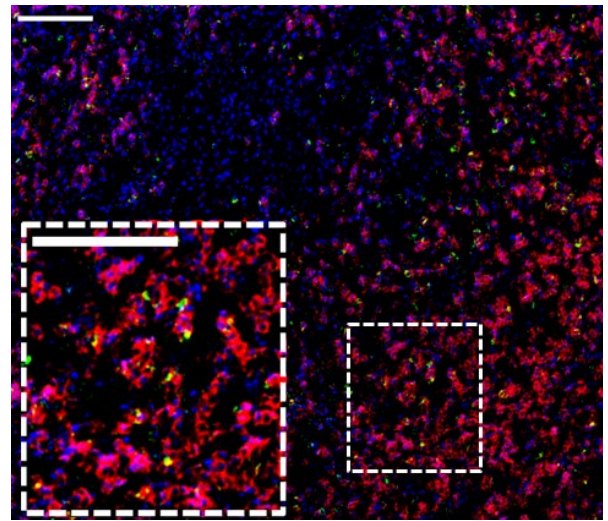

Human anaplastic large cell lymphoma (FFPE) stained with 153Eu-anti-CD223 (D2G40) at a dilution of 1:50 (green pseudocolor), 170Er-anti-CD3 (poly) (red pseudocolor), and iridium DNA intercalator (blue pseudocolor). Heat-mediated antigen retrieval was performed using Tris/EDTA buffer pH 9. Scale bar size = 100 µm.

## References

Chang, Q. et al. "Staining of frozen and formalin-fixed, paraffin-embedded tissues with metal-labeled antibodies for imaging mass cytometry analysis." *Current Protocols in Cytometry* 82 (2017): 12.47.1–12.47.8.

Giesen, C. et al. "Highly multiplexed imaging of tumor tissues with subcellular resolution by mass cytometry." *Nature Methods* 11 (2014): 417–22.

For technical support visit <http://techsupport.fluidigm.com>. | For general support visit [www.fluidigm.com/support](http://www.fluidigm.com/support).

### For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from CST™ and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information:** [www.fluidigm.com/sds](http://www.fluidigm.com/sds). **Patent and license information:** [www.fluidigm.com/legalnotices](http://www.fluidigm.com/legalnotices).

**Limited Use Label License:** The purchase of this product conveys to the purchaser the limited, non-transferable right to use the purchased consumable or reagent only with Fluidigm Instruments and Systems. **EU's WEEE directive information:** [www.fluidigm.com/compliance](http://www.fluidigm.com/compliance). Fluidigm, the Fluidigm logo, Hyperion, Imaging Mass Cytometry, and IMC are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. All other trademarks are the sole property of their respective owners. © 2017 Fluidigm Corporation. All rights reserved. 10/2017