

## Anti-Human IL-17A-172Yb

Catalog: 3172020B

Clone: BL168

Package Size: 100 tests

Isotype: IgG1

Storage: Store product at 4°C. Do not freeze.


Formulation: Antibody stabilizer with 0.05% Sodium Azide

Reactivity: Human

## Technical Information

**Validation:** Each lot of conjugated antibody is quality control tested by CyTOF® analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

**Recommended Usage:** The suggested use is 1 µl for up to 3 X 10<sup>6</sup> live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

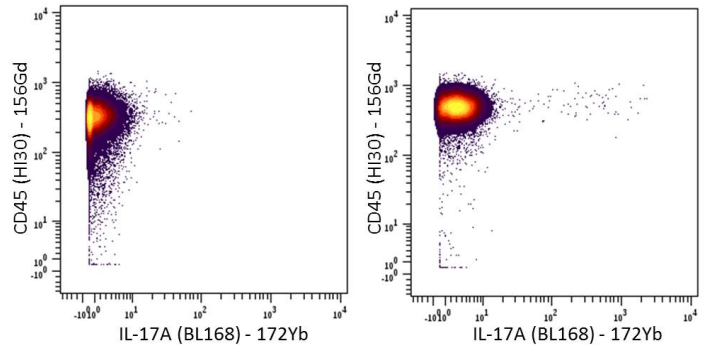

Human PBMCs were incubated for 6 hours in media alone (left) or with PMA and Ionomycin (right) in the presence of monensin and brefeldin A. Cells were then fixed, permeabilized, and stained with 156Gd-anti-CD45 (HI30) and 172Yb-anti-IL-17A (BL168).

## Description

Interleukin-17A (IL-17A) is a cytokine in a family that includes six structurally related proteins (IL-17A through IL-17F) though most activities seem to be mediated primarily by IL-17A. IL-17 induces neutrophil-rich inflammatory reactions, such as stimulating the production of chemokines and other cytokines, and stimulates the production of antimicrobial substances. It is primarily produced by a subset of T cells when activated (TH17 cells).

## References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. **Highly Multiparametric Analysis by Mass Cytometry.** *J Immunol Methods* 361 (1-2):1-20, 2010.

### For technical support visit [fluidigm.com/support](http://fluidigm.com/support)

EMAIL **United States and other countries not in Europe or Asia** [techsupport@fluidigm.com](mailto:techsupport@fluidigm.com) **Europe** [techsupporteurope@fluidigm.com](mailto:techsupporteurope@fluidigm.com) | **Asia** [techsupportasia@fluidigm.com](mailto:techsupportasia@fluidigm.com)  
 PHONE **United States (toll-free)** +1 866 358 4354 | **Europe** +33 1 60 92 42 40 | **Japan** +81 3 3662 2150 | **China** +86 21 3255 8368 | **All other countries** +1 650 266 6100

### For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information** [fluidigm.com/sds](http://fluidigm.com/sds) **Patent and license information** [fluidigm.com/legalnotices](http://fluidigm.com/legalnotices) | Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2015 Fluidigm Corporation. All rights reserved. 01/2015