

## Anti-Human IL-2-158Gd

**Catalog #:** 3158007B

**Package Size:** 100 tests

**Storage:** Store product at 4°C. Do not freeze.

**Cross Reactivity:** Human

**Clone:** MQ1-17H12


**Isotype:** Rat IgG2a

**Formulation:** Antibody stabilizer with 0.05% Sodium Azide

### Technical Information

**Validation:** Each lot of conjugated antibody is quality control tested by CyTOF® analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

**Recommended Usage:** The suggested use is 1 µl for up to 3 X 10<sup>6</sup> live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

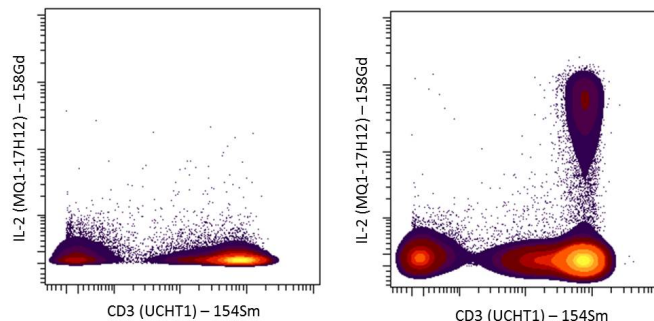

Human PBMCs were incubated for 5 hours in media alone (left) or with PMA and ionomycin (right) in the presence of monensin and brefeldin A. Cells were then fixed, permeabilized, and stained with 154Sm-anti-CD3 (UCHT1) and 158Gd-anti-IL-2 (MQ1-17H12).

### Description

Interleukin-2 (IL-2), originally known as T cell growth factor (TCGF), is a cytokine that plays a critical role in promoting and controlling T cell responses and functions. IL-2 stimulates the survival, proliferation and differentiation of antigen-activated T cells and is required for the survival and function of regulatory T cells. Clone MQ1-17H12 antibody is a non-neutralizing antibody.

### References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Newell, E.W., et al. Cytometry by Time-of-Flight Shows Combinatorial Cytokine Expression and Virus-Specific Cell Niches within a Continuum of CD8+ T Cell Phenotypes. *Immunity* 36 January 2012: 142-152.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. *J Immunol Methods* 361 (1-2):1-20, 2010.

### Contact Information:

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com  
www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.  
This product contains antibodies manufactured by and sold under license from eBioscience (XX-7029) and licensees thereof.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

