

Product Information Sheet

Anti-Mouse MHC Class I-144Nd

 Catalog #: 3144016B
 Clone: 28-14-8

 Package Size: 100 tests
 Isotype: Mouse IgG2a


Storage: Store product at 4°C. Do not freeze. Formulation: Antibody stabilizer with 0.05% Sodium Azide

Cross Reactivity: Mouse

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 μ l for up to 3 X 10⁶ live cells in 100 μ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

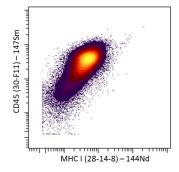

Mouse splenocytes stained with 147Sm-anti-CD45 (30-F11) and 144Nd-anti-MHC I (H-2Db) (28-14-8). Total viable cells are displayed in the analysis.

Description

The 28-14-8 monoclonal antibody reacts with the mouse MHC class I, H-2Db, and cross-reacts with H-2Ld, H-2Dq and/or H-2Lq. Binding of 28-14-8 to the alpha 3 domain of H-2Ld is not dependent on β2-microglobulin.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. Analytical Chemistry 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. J Immunol Methods 361 (1-2):1-20, 2010.

Contact Information:

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.

This product contains antibodies supplied by eBioscience.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

